

SPECIALE TITEL HER!

Anders Holst Rasmussen

May 18, 2021

1 Introduction

Here is a nice introduction.

2 Experimental Methods

The data used in this thesis stems from an experiment carried out at the IGISOL facility, at the University of Jyväskylä. Originally the experiment was scheduled for May 2020, but was postponed due to the ongoing Corona pandemic.

2.1 Detector setup

The detection setup consists of 6 double sided silicon detectors (DSSD), and 6 single sided silicon detectors (SSD). The detectors are $5\text{ cm} \times 5\text{ cm}$ and placed in a cube around the a target in the center, as shown on figure ???. The setup is designed with measuring the opening angle of the β in mind, and therefore the detectors covers 51% of the solid angle.

2.2 Experimental setup

The setup is designed to measure $\beta\alpha$ angular correlations in the β -delayed particle decay of ^8Li . When measuring multiple particles, the setup is highly dependent on the coverage of the solid angle. Therefore the setup is designed to have a large solid angle coverage, with high α -particle resolution, while still being able to measure β -particles.

This is has been achieved by creating a cube of six double sided silicon detectors (DSSD), all backed by a 1 mm single sided detector (SSD). To gain the largest solid angle, the detectors where placed as close to one another as possible. A 3D printed case was designed to hold the detectors in place, and achieved a solid angle coverage of 51% for the DSSD's. An illustration of the setup, together with the different detectors' thickness can be seen on

fig. 2.1. Even though the setup was designed to hold 12 detectors in total, there were only 11 detectors in the actual experiment. This was due to one of the SSD's being defect, so it was removed from the setup.

2.3 The detectors

As mentioned above, there were two types of detectors present in the setup. The first type is the Double sided silicon detector. As the name suggests, it consists of two sides, a front layer and a back layer. Each layer consists of 16 strips, that are placed in rows next to each other. The two layers are then arranged so each side are mutually orthogonal, which defectively makes pixels where each strip intersects a strip on the other side.

The strips on the front side are p-doped, while the back side are n-doped. When a charged particle hits the detector, it will ionize the atoms in the semi-conductor, and produce an electron-hole pair. The number of electron-hole pairs is proportional to the energy of the charged particle. The bias voltage on the detector collects the electrons and holes on opposite sides of the strip, where the charge is collected on aluminum contacts and a signal is measured. Energy is not deposited in these contacts, and therefore they constitute to a so called dead layer.

The detectors are square 5×5 cm and with their 16×16 strips, they have an effective grid of 256 pixels of 9 mm. 4 of the 6 detectors have a thickness of $60 \mu\text{m}$ and a dead layer of 100 nm dead layer. These detectors are the ones called Det1, Det3, Det4 and DetU in the setup seen on fig. 2.1. The other 2 detectors (Det2 and DetD) were both 1 mm.

The other type of detector was the SSD's also known as pads. They are different from the DSSD's, in that they only have one side, and no strips. Therefore they do not contain a grid, the same way the DSSD's do, and will therefore not provide any information as to where on the detector a particle has hit. But the reason they are in this setup, is because they are thicker. They are used to measure the energy deposited by a β -particle, and to determine if a given particle is a β -particle. An α -particle will deposit all of

its energy into the first DSSD it encounter, and be effectively stopped completely by it. But a β -particle will deposit almost no energy in a thin DSSD, and travel through it, ending up depositing some energy in the SSD. Therefore one can roughly distinguish the α -particles from the β -particles, by observing whether a particle has a hit in the DSSD and the SSD.

Since there are also two thick DSSD's in the setup, one can also get some information from a β -particle from these thicker detectors, as it will deposit more energy in these detectors.

2.4 AUSAlib and ROOT

ROOT is an object oriented C++ framework that is designed primarily for data analysis in high-energy and nuclear physics. It was created at CERN in 1995, and has since grown and become the dominant analysis software at both CERN and many other nuclear and particle physics laboratories. ROOT was designed to handle large amounts of data with high computing efficiency.

ROOT makes an intelligent data structure by creating a "Tree" with the class `TTree`. This tree will then have "branches" which corresponds to some variable of the given detection event, such as the energy of the front strip or identity of the detector. This `TTree` then allows for reading of an individual branch, while ROOT takes care of the memory management. One can also store a `TTree` to the disk in the form a `.root` file.

AUSAlib is a tool that build on top of ROOT. It was created by the subatomic group at Aarhus University. Before this tool was created, everyone in the group had to more or less create their own tools to get data from the detectors into a useful data structure. This meant that a lot of time was wasted just trying to access data from experiments. AUSAlib was therefore created, so the basic tasks of data extraction was automated.

AUSAlib has a lot of functionalities, but there are **XX** tools that was used in this thesis, namely the *Sorter*, *Calibrator* and the *AbstractSortedAnalyzer*.

2.4.1 Unpacker

The **Unpacker** converts raw data from the detectors into a ROOT TTree. This is done by using the unpacking program **ucesb** [ref.](#) This will setup the branch structure of the data. Some of these branches are **FT** and **BT**, which is a vector of the TDC (time) values for each event, for the front and backside of the detector. They are vectors because they contain information for each particle hits in a given time slot, for which there can be multiple. There are also the branches **FE** and **BE**, which is the ADC (energy) associated with the events.

2.4.2 Calibrator

Since the detectors work, by measuring an electrical charge that comes from the charged particle, the detectors needs to translate a specific charge to energy deposited.

To do that, we use the **Calibrator**-tool, which is designed to convert a channel number into an actual energy. Assuming that the channel numbers are linearly related to the energies, a known radioactive source can be measured, and the expected spectrum can be compared to the measured. This is done for each strip in each detector.

The **Calibrator** starts by running a peak-finding algorithm over some calibration data, to roughly identify the locations of the peaks, followed by a multi-Gaussian fit to find the most precise peak location. The positions of the peaks can then be compared to the expected energies, giving an associated energy to a given channel.

As mentioned earlier, all of the detectors have a small aluminum dead layer. All particles that pass through this layer will loses some amount of energy depending on the stopping power of the material and the effective thickness of the dead layer Δx_{eff} , which furthermore depends on the angle of incidence, θ . The relationship between the effective thickness and the actual thickness

is described as $\Delta x_{eff} = \Delta x / \cos(\theta)$. This gives the measured energy as

$$E' = E - \frac{dE}{dx} \frac{\Delta x}{\cos(\theta)},$$

where E is the original energy of the particle, dE/dx is the stopping power of the material and Δx is the thickness of the dead layer.

These calculations are all handled by the **Calibrator**. As input it takes an unpacked measurement of a source, a file specifying the locations of the expected peaks and a file specifying the spacial locations of the detectors. From this it calculates the energy loss, and creates a linear relationship between channel numbers and energies. This is then written to the disk as a separate calibration file, which can be parsed to other modules.

It is important to note that the Calibrator does not modify any data. Therefore the energy loss is unaccounted for. Instead it corrects the expected energy spectrum, which means that the resulting calibration is still valid. The energy loss correction is therefore still needed in the analysis, as the effect is unaccounted for in measurements.

2.4.3 Sorter

The sorter is used after a successful calibration. It generates a ROOT file based on the unpacked data, and applies the calibration. It is also responsible for matching and combining events from the front-side to the back-side of the detector. If there were one hit in the front side and one in the back, the matching is fairly trivial. If there however were multiple hits in both front and back, the **Sorter** will run a matching algorithm, which pairs the hits with the lowest energy differences.

When the events have been matched, the hits on the individual sides of each detector are merged into a single event. Therefore each event can now be considered a multiple of particle hits. This makes it possible to associate physical properties with each particle, such as direction and energy.

There has still not been done any filtering of the data, which is what we will discuss in chapter [3](#).

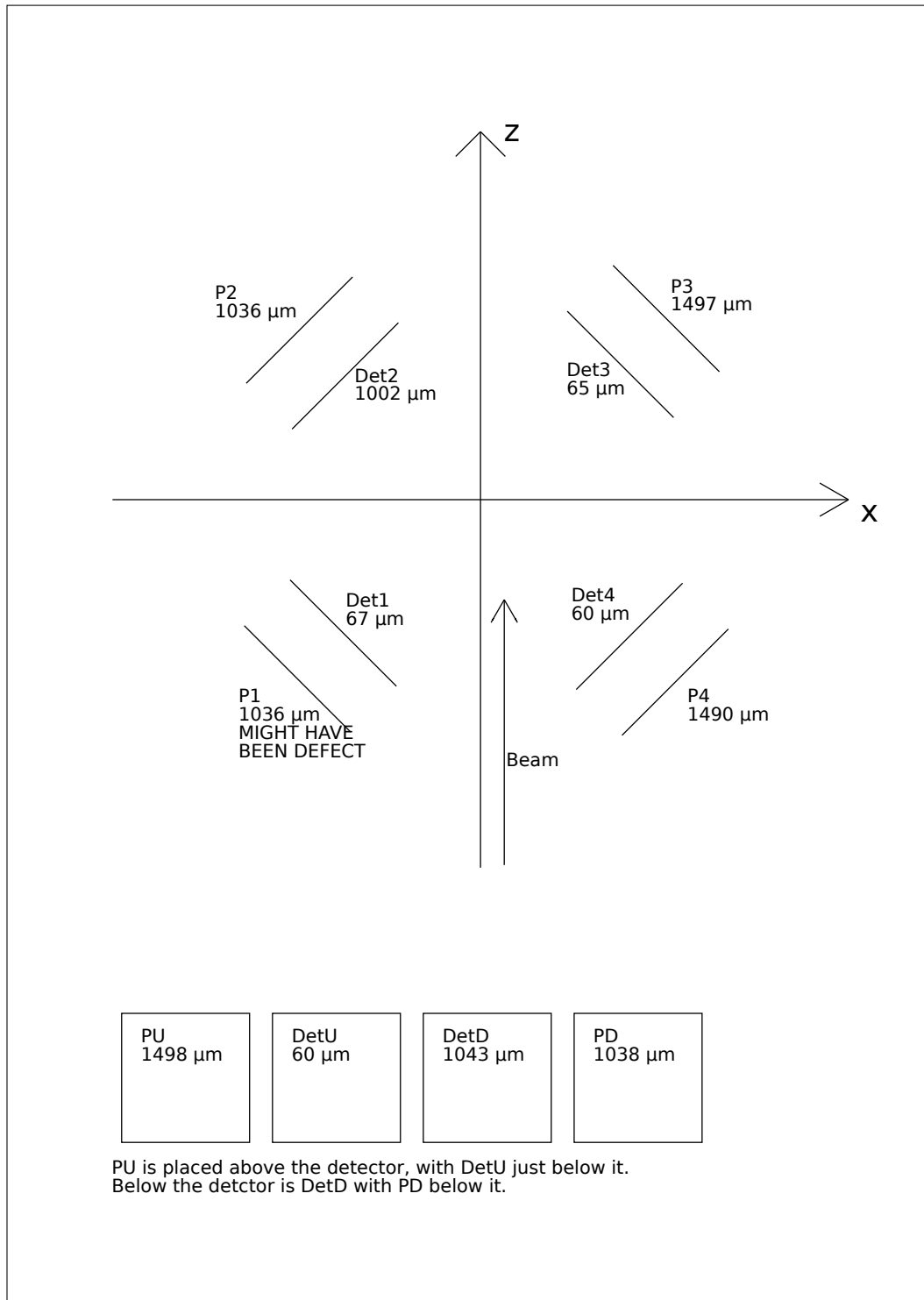


Figure 2.1: illustration of the setup. Note that the SSD P1 was absent from the actual experiment, due to failure in the detector

3 Data Reduction

3.1 Calibration

A typical single strip spectrum is shown on figure 3.1, where the calibrator has given an estimate of where the peaks are, illustrated by the red triangles. The radioactive source used to calibrate this setup contained ^{148}Gd , ^{239}Pu and ^{244}Cm . The proprieties of which is listed in Table 3.1.

Isotope	E_α [keV]
^{148}Gd	3182.690
^{239}Pu	5105.5
	5144.3
	5156.59
^{244}Cm	5762.64
	5804.96

Table 3.1: Decay energies for each isotope used in the calibration.

To compare the spectra, we need to find the peaks of the measured spectrum.

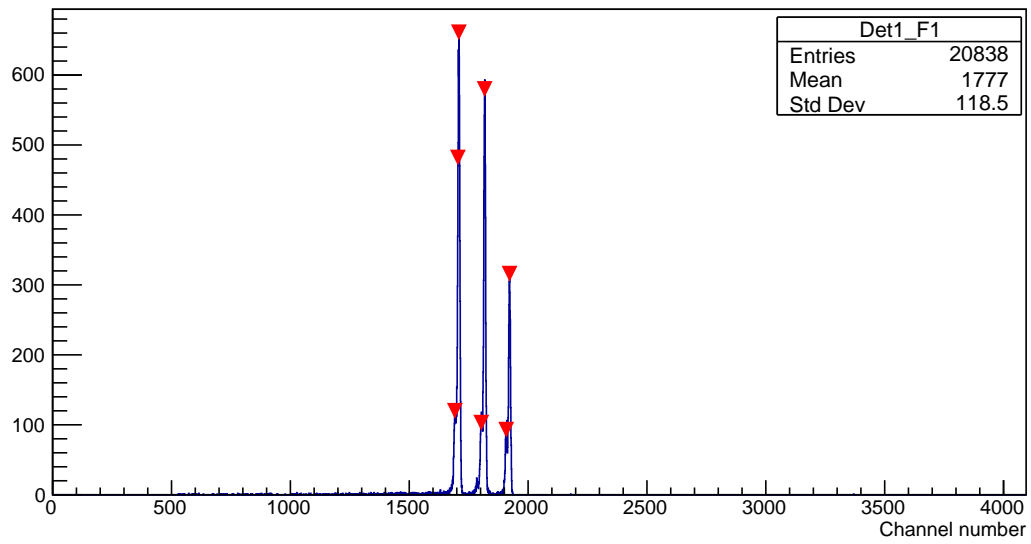


Figure 3.1: Calibration spectra from a front strip in Det1.